

Perforin Antibody : B-D48 PE conjugated

ANTI HUMAN PERFORIN PE CONJUGATED

CATALOGUE N° : 854.952.010

CLONE : B-D48

ISOTYPE : Mouse IgG1

IMMUNISATION : Natural human perforin

HYBRIDOMA : Myeloma X63/AG.8653 x Balb/c spleen cells

SPECIFICITY : Recognizes the human perforin, a 70 kDa protein

QUANTITY : 100 tests / 1ml


PURIFICATION : Ion exchange chromatography

FORMAT : Phosphate-buffered saline with 5% BSA and 0.1% sodium azide


WORKING DILUTION : Reconstitute with 1 ml deionised water. Use 10 µl to label 10⁶ cells or 100 µl of whole blood

STORAGE : Stable at +2-8°C for six months after reconstitution. DO NOT FREEZE

APPLICATION : Suitable for flow cytometry

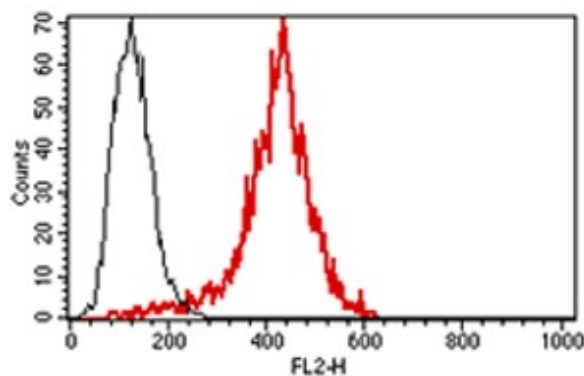

Flow cytometric analysis of perforin in YT-Indy cell line with the B-D48 monoclonal antibody

 Sodium azide is a poisonous and hazardous substance, handle and dispose of properly